


Rotational maintenance system guarantees thorough cleaning and lubrication of each handpiece. Coolant water and air channels are flushed.

### Low Speed Handpieces

Low speed handpieces incorporate gearing systems which require extra solution penetration. The NSK rotation system guarantees optimum oil delivery to the internal components of all low speed handpieces, while particles in the coolant water and air channels are removed by the cleaning solution.

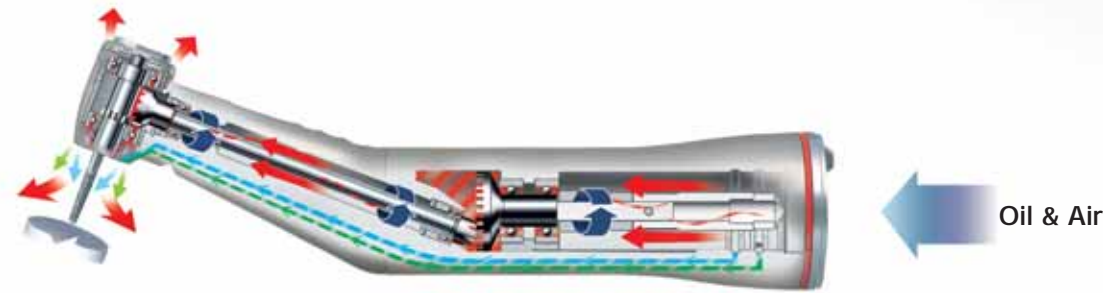

### High Speed Handpieces

The critical turbine component, plus all air and water channels of high speed handpieces are thoroughly cleaned and lubricated.




**SELECT** Two types of units are available according to individual demand.

**C1 Type** Single Rotation Port + Dual Standard Ports



**C2 Type** Dual Rotation Ports + Single Standard Port

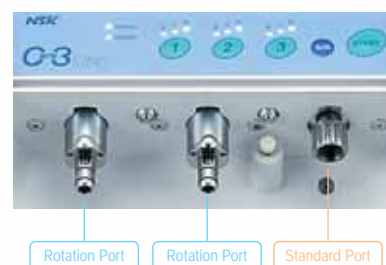

### Care3 Plus Complete Pack



#### Single Rotation Port + Dual Standard Ports

MODEL	ORDER CODE
Care3 Plus C1 Type (120V)	Y141-202
Care3 Plus C1 Type (230V)	Y141-204

#### Dual Rotation Ports + Single Standard Port

Care3 Plus C2 Type (120V)	Y141-203
Care3 Plus C2 Type (230V)	Y141-205

#### Contents

- Care3 Plus unit
- Power cord with L-shaped plug
- Maintenance oil (1 liter)
- Oil replenish funnel with filter
- Mist filter set
- Air tube
- Air filter

#### Specifications

- Power Supply : AC 120V or 230V 50/60Hz
- Air Pressure : 0.35-0.6 MPa (3.5-6 kgf/cm<sup>2</sup>)
- Weight : 7.1 kg (C1 type), 7.2 kg (C2 type)
- Dimension : W280 x D275 x H360 mm
- Solution Tank Capacity : 350 mL

- ◆ In full tank about 3,500 handpieces (S mode) are cleaned and lubricated.
- \* Connects only to Midwest 4-hole type
- \* Air piping is necessary (3.5-6 kgf/cm<sup>2</sup>)

### Care3 Plus Adaptors

MODEL	Care3 Plus ADAPTORS	ORDER CODE	MODEL	Care3 Plus ADAPTORS	ORDER CODE
E-Type Adaptor	For ISO E-Type & E-Type LUX Attachments	U433-035	SR Adaptor	For Sirona® High Speed Handpieces for Quick Coupling	Z257-023
PTL Adaptor	For High Speed Handpieces for PTL-CL-4HV-T Coupling	Z257-010	WH Adaptor	For W&H® High Speed Handpieces for Roto Quick® Coupling	Z257-021
QD Adaptor	For MACH-QD & CH-QD High Speed Handpieces	Z257-011	BA Adaptor*	For Bien-Air® High Speed Handpieces for Unifix® Coupling	Z280-022
B2/M4 Adaptor	For Borden 2-hole High Speed Handpieces	T904	F-Type Adaptor*	For ENDO-MATE TC, DT/TASKAL Wizard Heads	Z280-042
KV Adaptor	For KaVo® High Speed Handpieces for MULTiflex LUX® Coupling	Z257-020			

\*For Care3 Plus only. Don't use when connected to Care3. Check with the dealer for the adaptors to mount other manufactures' handpieces.

\*KaVo® and MULTiflex LUX® are registered trademarks of Kaltenbach & Voigt GmbH & Co., Germany. Sirona® is a registered trademark of Sirona Dental Systems GmbH, Germany. W&H® and Roto Quick® are registered trademarks of W&H Dentalwerk Burmoos GmbH, Austria. Bien-Air® and Unifix® are registered trademarks of Bien-Air Dental S.A., Switzerland.

### Spare Parts

MODEL	PRODUCTS	ORDER CODE	MODEL	PRODUCTS	ORDER CODE
Maintenance Oil (1 Bottle)	(1 liter)	Z016-117	AirTube (3m)		U508-038
Maintenance Oil Set (Pack of 6 Bottles)		Y900-138	Air Filter		U208-352
Mist Filter Set (pack of 12 pieces)		Y900-159	Air Filter Element		U508-353
Oil Replenish Funnel with Filter		U433-090	Fuse (120V)		U433-222
			Fuse (230V)		U433-223

**NAKANISHI INC.**  
www.nsk-inc.com

700 Shimohinata Kanuma-shi  
Tochigi 322-8666,  
Japan

TEL : +81 (0)289-64-3380  
FAX : +81 (0)289-62-5636

**NSK Europe GmbH**  
www.nsk-europe.de

Westerbachstraße 58  
D-60489 Frankfurt,  
Germany

TEL : +49 (0)69 74 22 99 0  
FAX : +49 (0)69 74 22 99 29

**NSK United Kingdom Ltd.**  
www.nsk-inc.com

Gateway1000, Arlington Business Park,  
White Way, Slough,  
Hertfordshire, SG1 2FP, UK

TEL : +44 800 634 1909  
FAX : +44 800 634 1910

**NSK Oceania Pty Ltd.**  
www.nsk-inc.com

Unit 22, 198-222 Young St.  
Waterloo, Sydney,  
NSW 2017 Australia

TEL : +61 (0)2 8306 3000  
FAX : +61 (0)2 8306 3030

**NSK America - Canadian Branch**  
www.nsk-inc.com

1054 Centre St., Suite 283  
Thornhill, Ontario, Canada  
L4J 8E5

TEL/FAX : +1 905 780 6757  
Toll Free : +1 866 780 6757

**NSK Middle East**  
www.nsk-inc.com

Dubai Airport Free Zone,  
Room 21W-310,  
3F West Wing, Dubai, UAE

TEL : +971 (0)4-299-6644  
FAX : +971 (0)4-299-6664

Specifications are subject to change without notice.

PR-D293E Ver.3 07.03.05.Y

**NSK**

AUTOMATIC HANDPIECE CLEANING & LUBRICATION SYSTEM

Care3 Plus

correct handpiece maintenance - everytime



# Care3 Plus Automatic Care3 Plus Is the Key to Handpiece Durability

Fast, user friendly, fully automatic, and economical.



Correct maintenance performed at regular intervals is a key to high performance and longer life of your handpieces. NSK designed easy-to-use Care3 Plus to facilitate handpiece maintenance at your office. Care3 Plus automatically cleans and lubricates up to 3 handpieces at one time with a touch of button. Both high speed and low speed handpieces are processed thoroughly in a single 2-minute operation cycle.



## Air Activation

Air supply is activated by the "Air" key. This enables the operator to flush out all residual solution.



## Original Air Filter

Clean air supply is a must for effective handpiece maintenance. NSK Care3 Plus comes with a highly efficient air filter that can be exchanged and installed easily wherever you desire.



## The Oil Level Indicator

Oil level indicator is ideally located at the front panel of Care3 Plus unit to clearly show amount of oil in the tank.



## Easy to Clean

At the touch of a button, the front cover can be completely removed for easy cleaning.



## Virtually No Mist

Spray mist is generally absorbed within the Care3 Plus unit by disposable mist filtration system. Filter sponges are easily exchanged. (Comes in a pack of 12)




## Power Cord with L-Shaped Plug

The power cord with L-shaped plug allows the Care3 Plus to be installed in virtually any surgery location.




Care3 Plus control panel is easy to understand and simple to operate.

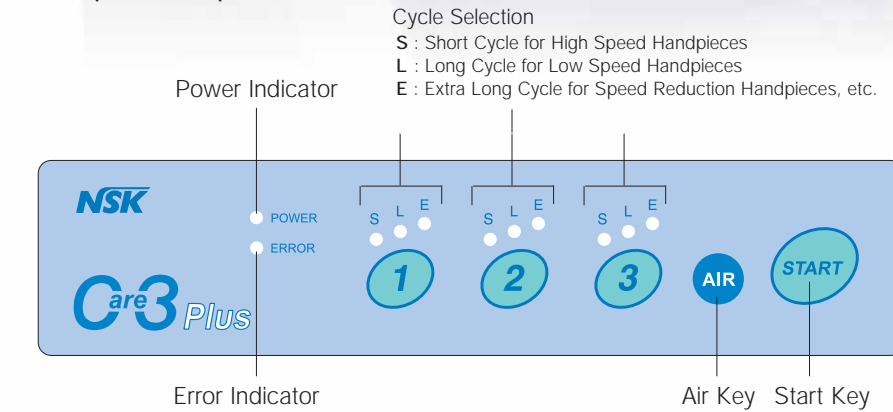

## Solution Tank

The solution tank has capacity for maintaining 3,500 handpieces (S Mode) before a refill is required. To refill the tank, simply unscrew the cap on top of the unit and pour in new Care3 Plus solution.

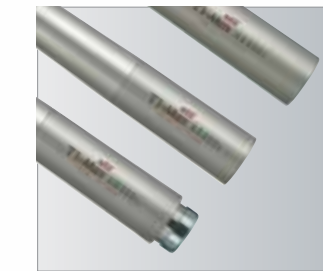

## Accommodates NSK and Other Brands' Handpieces

NSK and all other brands of E-Type low speed handpieces connect directly to the coupling supplied. NSK high speed handpieces require the PTL adaptor. Adaptors are available for direct connection of most other brands of high speed handpieces.

## Convenient Memory Function

Care3 Plus memory system stores and instantly recalls the last operational setting when the power switch is turned on. This memory function eliminates the need to repeatedly preset handpiece maintenance.

## End-of-Cycle Indicator

An audible indicator alerts the operator when the cycle is completed.

## Simple to Operate



1 Open the door and attach up to 3 handpieces to the Care3 Plus couplings.



2 After connecting the handpieces, close the door. Care3 Plus will not operate if the door is not properly closed.



3 Select the S, L, or E mode relevant to each handpiece connection location.

**S** : S mode (Approx. 35 seconds)  
For high speed handpieces

**L** : L mode (Approx. 45 seconds)  
For low speed handpieces

**E** : E mode (Approx. 50 seconds)  
For speed reduction and speed increasing handpieces, etc.

\*The figures in brackets indicate the clean-lubricate-flush time for a handpiece.



4 Press the start key. In less than two minutes, all handpieces are correctly cleaned and lubricated.

(When one long mode and two short modes are selected.)