

# Practical, compact and extremely reliable


X-RAY FILM DEVELOPER DÜRR XR 04



The wall-mounted solution for intraoral X-ray developing



*XR 04 develops intraoral X-ray films in a perfect image quality.*

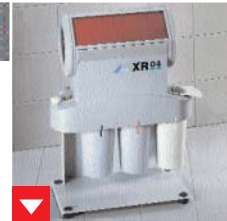

*XR 04 can be installed on the wall, on the table or in a cupboard.*



*Dür X-ray chemistry is specially adapted to the XR 04 guaranteeing best results.*

## The cost and time-saving alternative to manual developing

Intraoral X-ray images are an important part of an accurate diagnosis, but only if the image quality is right. Faults in film development make X-rays unusable. Dür Dental has therefore designed an X-ray film developing machine which is simple and time-saving to operate and produces excellent results, the XR 04 X-ray film developer. The first step on the road from manual developing to the automatic developing process. For consistently good developing results.

- ▶ Develops intraoral format 2 x 3, 3 x 4 and 2.7 x 5.4 cm X-ray films
- ▶ Easy to use, simple care and maintenance
- ▶ Includes integrated daylight attachment making the dark room unnecessary
- ▶ Variable processing time
- ▶ With or without heating

**DÜRR  
DENTAL**

# Save the valuable time of your assistants

## X-RAY FILM DEVELOPER DÜRR XR 04

### What could be more important to your assistants?

Time is an important factor in a dental practice: patients have to be looked after, there is paperwork to be done, the telephone rings – fine if one can count on a reliable X-ray developer, such as the XR 04 from Dürr Dental.



### Three steps to the perfect X-ray image.

Everything superfluous was deliberately dispensed with in the development of the XR 04, and therefore it is so safe and simple to operate. The exposed X-ray film is unpacked in the integrated dark room and loaded into the machine, the film then passes exactly through the various immersion baths, thanks to the perfected transporting system, and after no more than 6 minutes the X-rays are ready for you in the film outlet. For you this means no streaking, no clouding, but rather, optimally developed images perfect for diagnosis.



### Problem-free cleaning and care.

Dürr Dental's XR 04 is also simple and time-saving in terms of cleaning. The containers and transporting line



can be effortlessly removed in only a few steps, which for you means saving time and costs, and actively taking the load off the practice team.

### It all depends on the chemistry.

Perfect image quality is always dependent on the composition of the chemicals. Dürr Dental's O4X-ray chemistry is specially adapted to the XR 04 and thus guarantees the best results. Developer and fixer are prepared directly in the appropriate beaker.



Technical Data	XR 04
1740-01 (without heating)	220 - 240 V, 50/60 Hz
1740-03 (with heating)	220 - 240 V, 50/60 Hz
Shortest processing time	1'50"
Heating	25°C
Heating up rate	1°C/3 min
Chemical durability	3 weeks or 200 films
Chemistry	O4X-ray chemistry
Dimensions D x H x B	210 x 540 x 360 mm
Weight	7.5 kg

DÜRR DENTAL AG  
Höpfigheimer Strasse 17  
D-74321 Bietigheim-Bissingen  
Germany

www.duerr.de  
info@duerr.de

