

AUTOCLAVE STERILINE⁺

PROFESSIONAL INTELLIGENT STERILIZATION





CONFIDENT
SALES INDIA PVT LTD



23-Liters



confidentonline.com

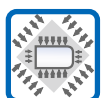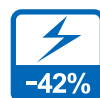



Normal working condition
Ambient temperature range : +5°C~+40°C
Relative humidity range : ≤85% ;
Atmospheric pressure : 70KPa~106KPa;
Water : Distilled water/Purified water.



German Ebro
Using German ebro temperature detector to check the temperature of each equipment, can meet the international European Union B class standard, temperature control +/-1°C, make sterilizing more thorough, more drier.



Solenoid Valve
Adopt Italian CEME original Solenoid Valve, large flow aperture, high temperature resistance, can bear pressure more than 6kgs, long service lift.



Door sealing self-adjustable
Bacterial Filter

Imported bacterial filter, filter air more efficiently while prolonging maintenance time.

Parameter

ORDER MODEL	MODEL	VOLTAGE	POWER	VOLUMN	CHAMBER SIZE	WEIGHT
MAX-8X	MAX-18L-III	110V-240V/50Hz-60Hz	1500W	18L	€245*355mm	49KG
MAX-8X	MAX-23L-III	110V-240V/50Hz-60Hz	2000W	23L	€245*465mm	54KG



LCD Display
LCD display function complete, can set up 15 national languages to meet customer needs, with power off after the prompt, machine fault prompt, and fault memory, convenient engineer home maintenance. Automatic atmospheric press detection adjustment, history query, unauthorized disabled etc.



Super Vacuum Pump
The working efficiency is fast, the limit vacuum is up to -90KPa, and the temperature resistance is up to 90°C.




Imported Vacuum Pump
Optional, Japan imported mute Vacuum Pump, light sound, circulating vacuum Principle, small Volume, fast working efficiency, stable.



Printer
In accordance with the regulation of the State Drug Administration, Print disinfection data function, reflect date, time, operation data(like working time, pressure, temperature etc.)



Over-current Circuit Breaker
Current unstable, automatic tripping, protect equipment

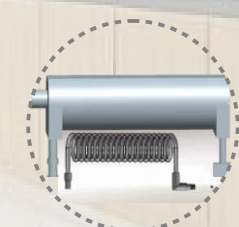

Innovative Steam Generator
The new steam generator concept combine with optimized program to solve the old steam generator blocking, seal aging and leakage problem. The resulting steam safe and reliable. Meantime the use of steam heating water tank, let the outlet water temperature increase to 45degrees or above. To speed up the sterilization time.(steam generator without cleaning all service life. German technology, Domestic leading)

Fast procedure

Min
(19-24)

wrapped instrument

(25-40)

Rubber material

(30-45)

Bloodthirsty procedure

(35-50)

Solid custom

(8-15)

Liquid custom

(15)

BD test

(25-40)

Vacuum test

(18-25)

HELIX test

HELIX test

(30-45)

Sterilization instrument

Unwrapped solid instrument, like pliers, holding needle, cleaning teeth and tweezers
This cycle is used for the metal instrument whose weight is below 1.0KG, like pliers, holding needle, cleaning teeth and tweezers

Fast procedure

Wrapped or unwrapped hollow instrument, Wrapped or unwrapped solid instrument
This cycle is used for solid and containing cavity metal instruments(weight:1.5kgs)with sterilizing high-pressure bag package (paper/paper, paper/plastic), also dental handpieces etc.

Rubber material

Rubber, fabric, wrapped hollow instruments, and so on.
This cycle is used for solid or unwrapped containing cavity metal instrument whose weight is 0.4kg.

Bloodthirsty program

hollow instruments which contains a bloodthirsty virus, solid instruments, dental cotton and so on.

Solid custom

Rubber, fabric, wrapped solid instruments, wrapped hollow instruments and so on.

Liquid custom

Glass and liquid

BD test

BD test pack, test the deposition uniformity of sterilizer steam temperature.

Vacuum test

Sterilizer sealing test

HELIX test

PCD challenge device, Steam sterilizer penetrability test



CONFIDENT
SALES INDIA PVT LTD

CONFIDENT SALES INDIA PVT. LTD.

No.47, 3rd Floor, Pete Chennappa Indl. Estate,
Magadi Main Road, Kamakshipalya, Bangalore – 560 079

Tel : 080-23002500 / 510, Fax: 080-23482156

E-mail : sales_enquiry@confidentonline.com

Toll Free : 1800 425 3366



confidentonline.com